

Lucira ankaptik COVID-19 Test Kit pitekn ta'n keknuwa'tikek

Kite'n msit ta'n ewikasik ke'smna'q poqtamkiaq

| Test kit ta'n koqowey pitek | Ke'smna'q poqtamkianuk |
|-----------------------------|---|
| | <p data-bbox="765 599 986 632">Waqma'ten kpitnn</p> <p data-bbox="1204 599 1486 632">Wi'ke'n wisun test- ituk</p> |

Lukwatis ta'n saqte'k. Kwile'n ta'n ma wen like aqq atestmulnin ta'n tujiw eskman 30 minutes test. Tlwula'tike aqq pipanikesi ta'n tujiw kitnman majukwatmn ta'n teluwek. Jiptuk ki'l tlo'tis l'nu'jaqmati wela'sik.

1 Ikanluwate'n wjit ilajuktmin test

| | | | | | |
|--|---|--|---|--|--|
| <p data-bbox="144 1240 338 1311">Panta'te'n packages 1 aqq 2</p> | <p data-bbox="421 1240 581 1311">Pija'te'n tapu'kl AA batteries</p> | <p data-bbox="690 984 878 1013">Amuj pa saqte'k</p> <p data-bbox="699 1212 852 1311">Ankapte'n stoqnamu'kn wasoqe'qn</p> | <p data-bbox="933 1240 1067 1274">Mna'te'n seal</p> | <p data-bbox="1159 984 1272 1013">menaqaj</p> <p data-bbox="1124 1240 1312 1274">Mukk paskija'tu</p> | <p data-bbox="1367 1212 1546 1311">Mukk panta'tu swap misoqo kuwayjk e'wmn</p> |
|--|---|--|---|--|--|



2 Sisqon ika'te'n swab

| | | |
|--|---|---|
| <p data-bbox="192 1654 378 1687">Mukk paskija'tu</p> | <p data-bbox="465 1404 633 1437">Mukk sama'tu</p> | <p data-bbox="1038 1393 1223 1581">1 2 Kasipute'n te's sisqon etu'k</p> <p data-bbox="847 1621 1442 1692">Msit swab pija'tu'n misoqo mist pitek ksisqonk Jiptuk tlo'tis wikawe'k.</p> |
|--|---|---|

3 Kito'qa'te'n swab vial


| | | | |
|---|--|---|---|
| <p data-bbox="126 1824 303 1891">Siawa'te'n ta'n metoqiaq</p> | <p data-bbox="581 1791 624 1824">15x</p> | <p data-bbox="847 1824 968 1891">Jikleketisk swab</p> | <p data-bbox="1310 2057 1571 2090">Siawa'sik piluwey page</p> |
|---|--|---|---|

4 Ankapte'n wula'sin test

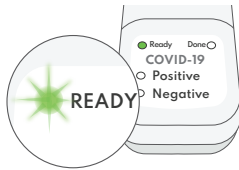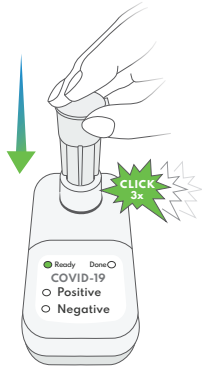

KEKNUWA'TASIK wasoqitesk te's 5 seconds. Katu mu wesogatesknuk, paskija'te'n ma'minu vial.

Skmate'n 30 minutes

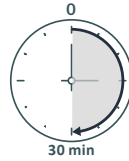

Pasik MU aja'tu test




Ta'n tujiw kaqkisa'tu'n, la'tu ta'n maliaptik airline wjit ta'n teluwek



Majukwate'n airline ta'n elukutijik telimiskik wjit ta'n teli jikleken test kit


Keknuwa'tikek

? Mu kelu'knuk kisa'tasik

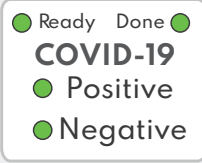


Pipanikesites piluwey kit aqq apa suwa'tis test

- Teluwek mu weskunumu'n

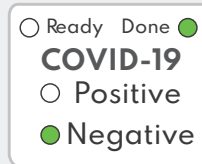

Nike' kisi siawa'sitis tujiw lietis ta'n pemilien Ula teluwek pasik eliklu'k 72 hours

+ Teluwek weskunmin

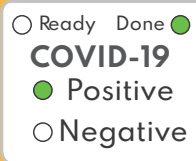

COVID-19 we'jitasik. Mukk we'kwata'siw, aqq knuwa'tu ta'n ikana'toq airline amkmayu. Jiptuk ma kisi liwun ta'n pemi lien.

Jiptuk amuj pa majukwatis wejiaq airport staff, airline staff, aqq ta'n pema'tu'tij wlo'ti wjit ta'n tlukwetis.

